



# PRESSURE CARE SOLUTIONS

FOR HOSPITAL ENVIRONMENTS



**Patient Comfort** 

Pressure Therapy

Staff Health & Safety







### **ACTIVE AIR® Bari 10** - Alternating and Turning Mattress System

A complete care solution combining pressure area therapy and patient handling assistance. Intelligent sensor system allows effective therapy for typical and bariatric patient needs.



Single - PTM639500 King Single - PTM639510 Super King Single - PTM639520



#### **Dual Alternation**

Patented variable alternation cycles

### Static Side Walls

Improves patient safety and manual handling efficiency

**Intelligent High Flow Pump** Automatically switches to low-flow rates when alternating to improve patient sleep

#### Multi-Stretch Vapour **Permeable Cover**

Helps to reduce shear forces whilst improving skin microclimate

## **Respiratory Therapy** / Pulsation Mode

Cycles between low pressure and firm inflation to promote oxygen saturation and improved perfusion during deep breathing exercises

#### **Integrated Turning Therapy** Allows patient rotation up to 40° for

improved pressure care and to assist staff with patient handling requirements



a single asset can meet the broadest spectrum of patient needs

**High Flow** Rapid CPR Valve for typical users (up to 160kg) and bariatric users (up to 450kg) ensuring

#### **RAPID PUMP INFLATION**

**AutoLock** 

Control panel automatically locks if not

touched for 30 seconds, preventing accidental

changes to modes and settings.

Minimises therapy response times and staff downtime





## **AUDIBLE / VISUAL ALARMS**

Integrated fault sensor and service indicator improves uptime and therapy outcomes



AutoProtect

Intelligently ramps up air flow to the sacral region as bed backrest is raised



(Clean with up to 10,000ppm

Sodium Hypochlorite solution)





30kg - 450kg





Breathable & Waterproof

## INTEGRATED TURNING THERAPY

- Lateral (left/right side focused) inflation and alternation is designed to redistribute pressure while also assisting with patient handling:
  - o Programmable time and turn direction to meet individual patient needs
  - Alternation therapy continues during turning mode due to 1 in 4 alternation
  - o Compliments patient repositioning programs
- Improves staff workplace health and safety



10 - 90min TURN CYCLE

## STAFF SAFETY

- Turning assistance provides a positive benefit for staff and patients by:
  - o Reducing the forces required by staff to manually reposition a patient
  - o Improving comfort for the patient during the repositioning process
- Facilitating more efficient daily processes such as:
  - Positioning slings and slide sheets
  - Patient repositioning, clothing, bathing and toileting
  - o Changing bed linen



raist-hip ratio, Report of a WHO expert consultation, Geneva, 8-11 December 2008

#### BARIATRIC BENEFITS


- Seated Fowler Position supports:
  - o Focused sacral therapy
  - Therapy continuation at meal times
  - Respiratory considerations
  - Pulsation therapy to compliment respiratory therapy programs
  - Constant Low Pressure for pain management




**TRANSVERSE ALTERNATION** 



CHECKER BOARD ALTERNATION

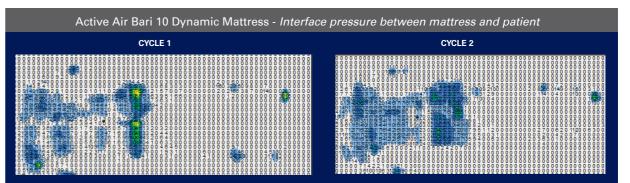

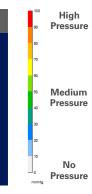



# ACTIVE AIR® Acute 8 with turnassist - Dynamic Mattress Replacement

The Active Air Acute 8 is specifically designed to meet the demands of Critical Care environments, maximising patient comfort and therapeutic outcomes.



Single - PTM639400T King Single - PTM639410T



RAPID INFLATOR COMPATIBLE





**AUDIBLE / VISUAL ALARMS** Integrated fault sensor and service indicator

9

turnassist valve



**Rapid Inflation Turbo Pump** 

(Optional accessory:

PTM639000)

**Hospital Grade** (Clean with up to 10,000ppm Sodium Hypochlorite solution)



Single | King Single





Breathable & Waterproof



Easy to operate valves located on the side of the mattress (foot end) engage lateral turning through integrated air cells. This function drastically reduces the WHS impact on care workers with a turning angle of up to 30° whilst also improving:

- In-bed washing
- Changing of Linen and Clothing
- Positioning of manual handling slings/sheets
- Targeted relief of asymmetric/laterally located pressure areas

## DESIGNED FOR ACUTE CARE ENVIRONMENTS



1 in 3 Cell Alternation Cycle with 24 Narrow Profile cells - improves patient comfort by reducing space between cells and increasing patient support surface area.




**Dual Compressor Pump** – reduces noise and efficiently manages patient therapy at highest patient weight loads (250kg).



Focused Heel Care - removable heel cells facilitate rapid response to changes in heel therapy needs.



Clear Communication – on pump patient weight guide, troubleshooting guides and mattress cover graphics provide critical information such as cleaning instructions, asset ID numbers and servicing considerations.

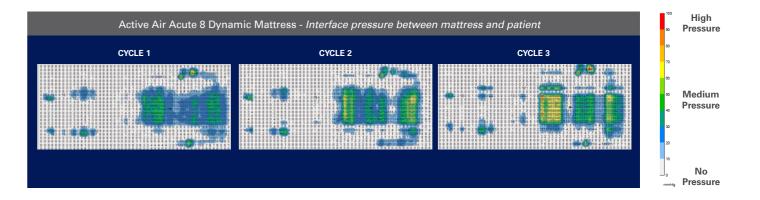




# **HYBRID** Comfort - Passive to Active Support Surfaces

The Hybrid Comfort is an Acute mattress system offering foam-in-cell technology enabling flexible operation as a high comfort foam static support surface or a powered-up alternating pressure support surface.



Single - PTM638005 King Single - PTM638006

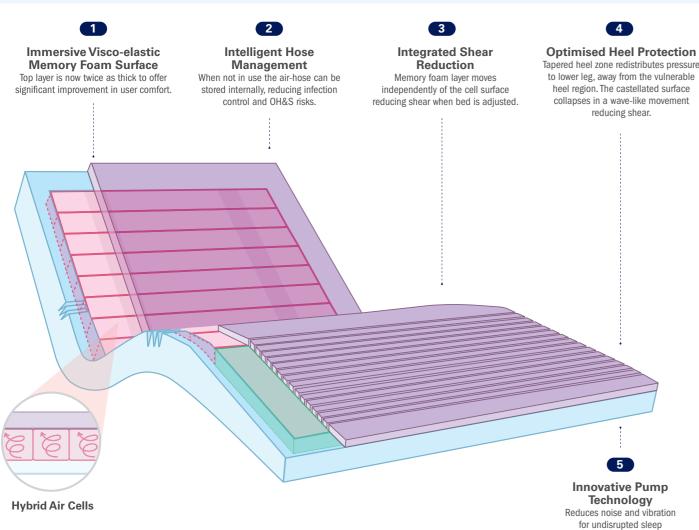

## **CLINICAL BENEFITS**

- Infection Control Fully welded seams, with waterfall zip covers, protect against fluid ingress
- Skin Microclimate Multi-stretch cover reduces shear and promotes air flow and optimal skin microclimate
- Sleep Quality Innovative pump technology reduces noise and vibration
- CPR Compatible In non-powered mode, the foam-core cells allow CPR to take place. When in powered mode, simply disconnect the pump to facilitate resuscitation
- Reduced Manual Handling Effort Mattress handles assist staff moving mattresses



**Hospital Grade** (Clean with up to 10,000ppm Sodium Hypochlorite solution)



Single | King Single








Waterproof

#### LIFECOMFORT®

## **Premier Care Acute**

**Engineered Multi-Layered Support** 

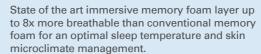

Single - BEM046735 King Single - BEM046740

Has three foam layers including a castellated patient surface providing users with immersive, pressure relieving support and air flow for enhanced skin microclimate

Castellated Static Support Surface


**Acute Care** 

LIFECOMFORT\*

Single - BEM046615 King Single - BEM046616



castellated cells, closely contours to the user's body for optimal pressure redistribution and comfort

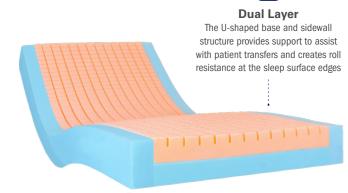

#### **Effective Manual Handling**

Integrated handles on both sides of the mattress surface assist safe manual handling when moving mattress between rooms and wards

and reduce manual handling burden on staff when assisting patients in and out of the bed

#### 4-Way Stretch Top Cover

4-way stretch top for shear mitigation and friction reduction and 4-sided zipper with waterfall covers that can be easily removed for effective cleaning and infection control

## Reinforced Side Walls

Maximise patient mobility

#### Clear Communication

Mattress cover graphics offer critical information such as cleaning instructions and space for asset ID numbers to be recorded

#### **Bed Surface Integration**

Hinged foam core allows mattress to move with profiling bed surfaces



Cleaning Compatible with solutions containing up to 10,000ppm Sodium Hypochlorite

#### **Mattress Rotation Guide**

Intuitive cover graphics remind nursing. caring and cleaning staff to continually rotate mattress for improved performance and longevity

## **CLINICAL BENEFITS**

- MEMair BREEZE™ top layer offers leading breathability for a cooler more immersive sleep surface
- Convoluted second layer coupled with varying foam grades

## **CLINICAL BENEFITS**

 Castellated patient surface promotes airflow through varying size channels. Sections move independently to redistribute shear and compression forces



**Hospital Grade** (Clean with up to 10,000ppm Sodium Hypochlorite solution)



Single | King Single



30kg - 250kg



Washable



**Breathable &** Waterproof





Active Air is a registered trademark.

# **SPECIFICATIONS**

	Active Air Bari 10	Active Air 8 Acute	Hybrid Comfort	Premier Care Acute	Acute Care
ACUTE ALTERNATIN	IG AIR & FOAM MATTRESSES				
Sizes & Codes	PTM639500 (Single) PTM639510 (King Single) PTM639520 (Super King Single)	PTM639400T (Single) PTM639410T (King Single)	PTM638005 (Single) PTM638006 (King Single)	BEM046735 (Single) BEM046740 (King Single)	BEM046615 (Single) BEM046616 (King Single)
System	Air Mattress Replacement	Air Mattress Replacement	Hybrid Mattress	Foam Mattress	Foam Mattress
Risk Category	High Risk - Very High Risk	High Risk - Very High Risk	Medium to High Risk	Medium to High Risk	Medium to High Risk
Patient Weight	30kg - 450kg	30kg - 250kg	50kg - 230kg	30kg-250kg	50kg-250kg
Number of Cells	20	24	9	N/A	N/A
Static Head Cells	0	3	N/A	N/A	N/A
Cell Cycle	1 in 2, 1 in 4	1 in 3	1 in 2	N/A	N/A
Cycle Time	10 / 15 / 20 mins	10 / 15 / 20 / 25 / 30 mins	11 mins	N/A	N/A
Cell Material	Pure TPU	Pure TPU	PU	N/A	N/A
DIMENSIONS (L x W x	H in mm)				
Single	2000 x 900 x 270	2000 x 900 x 200	1980 x 895 x 170	2000 x 900 x 150	1980 x 880 x 150
King Single	2000 x 1070 x 270	2000 x 1050 x 200	1980 x 1050 x 170	2000 x 1050 x 150	2030 x 1050 x 150
Super King Single	2000 x 1230 x 270	N/A	N/A	N/A	N/A
Weight	11kg - 16kg	9kg - 11kg	20kg - 22kg	14kg - 16kg	14kg - 16kg
Safe Working Load	450kg	250kg	230kg	250kg	250kg
ARTG Number	293518	293518	291660	288486	288486
PUMP					
Controls	Digital	Digital	Digital	N/A	N/A
Weight	6.5kg	4.6kg	3.5kg	N/A	N/A
Dimensions (mm)	400 x 175 x 275	420 x 130 x 230	320 x 130 x 230	N/A	N/A
Supply Voltage	230V, 50Hz	230V, 50Hz	230V, 50Hz	N/A	N/A
FEATURES					
AutoDetect	N/A	✓	N/A	N/A	N/A
AutoProtect	✓	✓	N/A	N/A	N/A
AutoLock	✓	$\checkmark$	N/A	N/A	N/A
CPR Release	✓	✓	N/A	N/A	N/A
Lateral Turning	✓	✓	N/A	N/A	N/A
WARRANTY					
Foam	N/A	N/A	5 Years	5 Years	5 Years
Cover	2 Years	2 Years	2 Years	4 Years	4 Years
Air Mattress	2 Years	2 Years	N/A	N/A	N/A
Pump	2 Years	2 Years	2 Years	N/A	N/A